


Indicator:



Methyl orange (yellow to red)

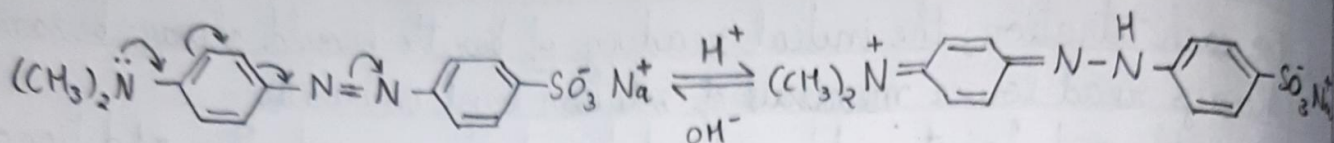



Benzenoid structure, yellow

Quinonoid structure, pink

Chemical reactions:

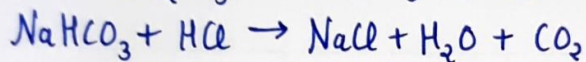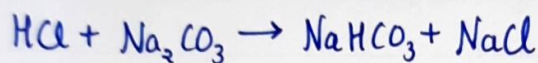

Ionic reactions:

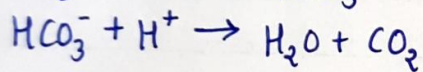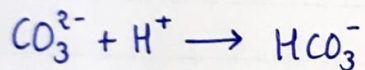

Observations:

Titration of given HCl v/s standard $\frac{N}{40}$ Na_2CO_3 solution (Table 3.1)

S. No.	Volume of Na_2CO_3 (ml)	Burette readings		Volume of HCl used (ml)	Concordant readings (V _o)
		Initial	Final		
1.	10 ml	0.0	6.0	6.0	6.0 ml
2.	10 ml	6.0	12.0	6.0	
3.	10 ml	12.0	18.0	6.0	
4.	10 ml	18.0	24.0	6.0	

Experiment 3

Aim:

To determine the total alkalinity of given water sample using methyl orange as indicator. Given standard $N/40$ Sodium Carbonate to standardize given HCl solution

Apparatus required:

Pipette, Burette, Stand, Conical flask, Measuring cylinder, White paper.

Reagents required:

Sodium Carbonate (Na_2CO_3), Hydrochloric Acid (HCl), Alkaline water sample, Distilled water.

Theory:

This is an example of neutralization reaction titration. The titration is based upon the specific pH range of the indicator. Methyl orange gives red colour in pH range 4.4-3.0. Initially when standard hydrochloric acid is added, CO_3^{2-} ions are converted to HCO_3^- . On addition of further acid, the HCO_3^- ions are converted to CO_2 and H_2O . When all HCO_3^- ions are used up the pH of solution becomes less than 4.4 instantaneously. At this point methyl orange gives red colour.

Procedure:

1. Standardization of given HCl solution
 - (i) Take 10 ml of standard Na_2CO_3 with the help of pipette into conical flask. Add 1-2 drops of methyl orange.



Titration of standardized HCl v/s given Alkaline water sample. (Table 3.2)

S.no.	Volume of water sample	Burette readings		Volume of HCl used (ml)	Concordant reading (V _o)
		Initial	Final		
1.	10 ml	0.0	10.5	10.5	10.5 ml
2.	10 ml	10.5	21	10.5	
3.	10 ml	21	31.5	10.5	
4.	10 ml	31.5	42.0	10.5	

Calculations:

To calculate Normality of given HCl solution, apply normality equation:

$$N_{HCl} \times V_{HCl} = N_{Na_2CO_3} \times V_{Na_2CO_3}$$

$$N_{HCl} = \frac{N_{Na_2CO_3} \times 10}{V_{HCl}} = \frac{10}{40} \times \frac{1}{6} = \frac{1}{24} N = 0.04167 N$$

To calculate Normality of given water solution, apply normality equation:

$$N_{sample} \times V_{sample} = N_{HCl} \times V_{HCl}$$

$$N_{sample} = \frac{N_{HCl} \times V_{HCl}}{10} = \frac{1}{24} \times \frac{10.5}{10} = 0.04375 N$$

Strength of Alkaline content in given water sample:

$$\begin{aligned} \text{Strength} &= \text{Normality} \times \text{Eq. wt. of } CaCO_3 \\ &= 0.04375 \times 50 = 2.1875 g L^{-1} \end{aligned}$$

$$\text{Strength} \approx 2.19 g L^{-1} \text{ or } 2187.5 ppm$$

- (i) Run the acid solution from the burette into conical flask drop wise with constant shaking of solution.
- (ii) Found the end point when light yellow colour solution turned red. Noted volume of acid used. Repeated this step 4-5 times till you got atleast 2 concordant readings.

2. Titration of given water sample with HCl solution

- (i) Pipette out 10 ml of given water sample into conical flask. Added 2-3 drops of methyl orange indicator.
- (ii) Added HCl solution from burette into the conical flask with constant shaking.
- (iii) Titrated with acid till yellow coloured solution turned red (end point). Noted the volume of acid used (V_2). Repeated the steps to get concordant reading.

Result:

Total alkalinity of the given solution (in terms of CaCO_3) is 2187.5 ppm

Precautions:

- (i) Before starting the experiment, the glass apparatus must be cleaned.
- (ii) For each titration, the initial readings of burette should always be same.
- (iii) Always read lower meniscus of solution level in burette.
- (iv) Near the end point, add acid solution dropwise and after addition of each drop, see the colour against white background.
- (v) Do not blow last drop of solution from pipette. Just tap the tip of pipette to the walls of the flask.

